- | | |
|---|---|
| (D) Hornhautentferner
Gebrauchsanweisung..... 2-7 | (RUS) Прибор для удаления роговицы
Инструкция по применению.....33-38 |
| (GB) Callus remover
Instructions for use..... 8-12 | (PL) Urządzenie do usuwania
zrogowaciałego naskórka
Instrukcja obsługi 39-43 |
| (F) Ponceuse
Mode d'emploi 13-17 | (NL) Eeltverwijderaar
Gebruiksaanwijzing 44-48 |
| (E) Eliminador de durezas
Instrucciones para el uso..... 18-22 | (P) Removedor de calosidades
Instruções de utilização 49-53 |
| (I) Levigatore per talloni
Istruzioni per l'uso 23-27 | (GR) Συσκευή αφαίρεσης κάλων
Οδηγίες χρήσης 54-60 |
| (TR) Nasır Sökücü
Kullanım kılavuzu 28-32 | |

Inhalt

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Contents

1. Intended use	9	5. Operation.....	12
2. Safety notes	9	6. Maintenance and storage.....	12
3. Unit description	11	7. Replacement parts and wearing parts...	12
4. Initial use	11		

Dear customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage, beauty, baby and air.

Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

With kind regards,
Your Beurer team

 WARNING

- The unit is only intended for domestic/private use, not for commercial use.
- This device may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- If the mains connection cable of this device is damaged, it must be replaced by the manufacturer, its customer service department or a similarly qualified person to avoid any risks.

Included in delivery

- Base device
- Coarse callus attachment
- Fine callus attachment
- Protective cap
- Protective cap for skin cells
- Charger
- Mains part plug
- Cleaning brush
- These instructions for use

Signs and symbols

The following symbols appear in the instructions for use and on the type plate:

WARNING Warning instruction indicating a risk of injury or damage to health.

IMPORTANT Safety note indicating possible damage to the unit/accessory.

Note Note on important information.

Note Suitable for cleaning under running water.

The device is double protected and therefore corresponds to protection class 2.

Only use indoors.

1. Intended use

With this callus remover you can remove unwanted hard or rough skin and calluses simply, quickly and effectively.

A coarse callus attachment and a fine callus attachment are included in the delivery, which you can use for optimum treatment of the affected areas.

The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

2. Safety notes

Read these instructions for use carefully. Non-observance of the following information may result in personal injury or material damage. Store these instructions for use and make them accessible to other users. Make sure you include these instructions for use when handing over the device to third parties.

WARNING

- Keep packaging material away from children! Choking hazard.
- Children must be supervised at all times whilst using the device.

IMPORTANT

- Before use, ensure that the device and accessories show no visible damage. If you have any doubts, do not use the device and contact your retailer or the specified Customer Services address.
- Do not pull, twist or bend the mains cable.
- Ensure that the mains cable is not routed or pulled over pointed or sharp objects.
- Insert the mains part plug into a socket that has been installed in compliance with applicable regulations.
- Never open the housing.
- Keep the device away from pointed or sharp objects.
- If the device has been dropped or has suffered any other damage, it must no longer be used.
- The device must be switched off and disconnected from the mains every time after use and before cleaning.

WARNING: Electric shock

As with every electrical device, use this device carefully and cautiously, in order to avoid hazards due to electric shock.

Therefore, observe the following instructions for use:

- Only operate the device with the mains part plug provided and the mains voltage specified on the mains part.
- Never use the device if it or its accessories show visible signs of damage.
- Do not use the device during pregnancy.

In the case of defects or operational faults, switch the device off immediately and disconnect the device from the mains socket. Do not pull on the mains cable or on the device to disconnect the mains part from the socket. Never hold or carry the device by the mains cable. Keep the cables away from hot surfaces.

Ensure that the charger, the plug and the cable do not come into contact with water or other liquids.

Therefore, observe the following instructions for use:

- Only use the charger in dry interior rooms.
- Only use the charger with dry hands.
- Never use the charger in the bath, shower or near a filled wash basin.
- Never use the charger in a swimming pool, jacuzzi or sauna.
- Never use the charger outdoors.

Never reach for a charger that has been submerged in water. Disconnect the mains plug immediately.

WARNING: Repairs

- Only specialist personnel may perform repairs on electrical devices. Improper repairs may subject users to considerable danger.
In the event of a fault or damage, have the device repaired by a certified specialist workshop.
- If the mains part or its cable is damaged, this must be replaced by an authorised retailer.

WARNING: Risk of fire

There is a risk of fire if the device is used improperly or these instructions for use are ignored.

Therefore, observe the following instructions for use:

- Always supervise the device when it is in use, especially if there are children nearby.
- Never use the device underneath a cover, such as blanket, pillow, etc.
- Never use the device near petrol or other highly flammable substances.

Note: Disposal


Observe the local regulations for material disposal. Please dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any queries, please contact the appropriate local authorities.

Remove the battery before disposing of the device at an official recycling point. Before removing the battery, disconnect the device from the charger and wait until the motor comes to a complete stop. Remove the front cover using a screwdriver and loosen the screws underneath. Finally, disconnect the battery. Dispose of the used battery at a collection point.

3. Unit description

1. Base device
2. Base station
3. Coarse callus attachment
4. Fine callus attachment
5. Unlock button
6. Fixing cap
7. LED light
8. Battery level indicator
9. On/off switch with 2 speed settings
10. Protective cap
11. Protective cap for skin cells

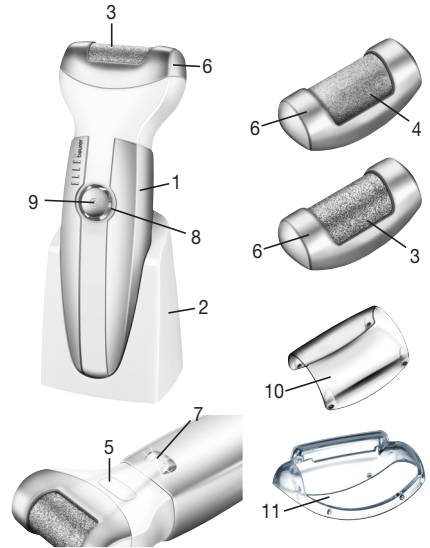4. Initial use

- Remove the packaging.
- Check the device, mains plug and cable for damage.
- Connect the device to the mains. When doing so, the device must be switched off.
- Ensure that the cable does not pose a trip hazard.
- Charge the device for 2 hours before the first use. If the battery level indicator flashes slowly this means that the battery is charging. Once the display is permanently illuminated, the device is fully charged. Then you can begin using the device either on battery power or using the mains cable.
- If the battery level indicator flashes rapidly during battery operation, the battery is empty. Charge the device by placing it on the charger or connecting it to the mains adapter.
- One charge is sufficient for approx. 30 minutes of cordless operation. After that you can continue to use the device via mains operation.

IMPORTANT

- For hygiene reasons, you should not share the attachments with anyone else.
- When switched on, the device should not come into contact with other areas of skin (e.g. eyelashes, head hair, etc.), or with clothes or cords, to avoid injury as well as blockages or damage to the device.
- Do not use the device on irritated skin, varicose veins, rashes, spots, (hairy) birthmarks or wounds without consulting a doctor beforehand. The same applies to patients with weakened immune systems, diabetes, haemophilia and immune insufficiency.

Note

- Batteries do not have their full load capacity from the beginning. In order to achieve the full load capacity of the batteries, they should be fully charged, completely discharged and then fully charged once more prior to initial operation of the device. This recharge/discharge cycle should be repeated as often as possible to extend the service life of the batteries.
- The callus remover can be used on dry or wet skin. Removal of calluses is more effective on dry skin.
- When used on dry skin, the skin must be completely dry and free from grease.
- Do not use any cream before removing calluses.
- Do not remove too much of a callus as it may otherwise grow back again more strongly.

5. Operation

- Press the on/off switch to switch on the device.
- If you press the on/off switch again, the speed setting increases. The LED light ensures optimal lighting conditions so that you can see and remove calluses better.
- Move the device slowly and evenly over the callus. Make sure you hold the device vertically at a 90° angle to the callus. To achieve optimum results, it may help to move the device over the callus in different directions.

Changing the attachment:

- To change an attachment, press the unlock button and take the attachment off the device. You can remove the fixing cap from the callus attachment by pressing the attachment hooks outwards.
- If you wish to remove highly callused areas, use the coarse callus attachment.
- If you wish to remove slightly callused areas or to smooth your skin, use the fine callus attachment.
- Place the protective cap for skin cells on the selected callus attachment to safely collect skin cells.

6. Maintenance and storage

WARNING

- Unplug the device each time before cleaning.
- Clean the device using only the methods specified. Under no circumstances may liquid enter the device or the accessory.
- Clean the callus attachment with the cleaning brush. We recommend that you moisten this beforehand with alcohol to guarantee a hygienic clean. Brush the attachments out from behind while you turn the roller by hand.
- Make sure that the relevant callus attachment is completely dry before you put it back on the device.
- Do not use the device again until it is completely dry.
- Do not clean the device in the dishwasher.
- Do not use any abrasive cleaning products or hard brushes.

Clean the unplugged device once it has cooled down using a soft, slightly damp cloth and a mild cleaning agent. Then wipe dry with a soft, lint-free cloth.

7. Replacement parts and wearing parts

Replacement parts and wearing parts are available from the corresponding service address (according to the service address list). Please state the corresponding order number.

Replacement part	Item number and/or order number
Fine callus attachment with fixing cap	163.280
Coarse callus attachment with fixing cap	163.281

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